

The materials and information provided via this brochure are distributed internationally, but not all the products or services found on this brochure are necessarily available in your country or area. Please contact local office or agents to make sure of the availability of the relevant products or services.

Mindray Medical Netherlands B.V.
Drs.Wan Royenstraat 8, 3871 AN Hoevelaken, The Netherlands
Voice: (31-33) 254-4911
Fax: (31-33) 253-4280

Mindray Medical Italy S.r.l.
Via Leonardo DA Vinci, 158, 20090 Trezzano Sul Naviglio (MI), Italy
Voice: (39-02) 5737-401
Fax: (39-02) 5523-2018

Mindray Medical España S.L.
Avenida Manoteras 38, Edificio B, Bajo B008, 28050-Madrid Spain
Voice: (34-91) 392-3754
Fax: (34-91) 392-3756

Mindray Medical Germany GmbH
Zwischen den Bächen 4, 64625 Bensheim, Germany
Voice: (49) 6251-17524-0
Fax: (49) 6251-17524-20

Mindray (UK) Limited
Mindray House, Kingfisher Way, Hinchingsbrooke Business Park,
Huntingdon, Cambs, PE29 6FN
Voice: (44) 1480-416-840
Fax: (44) 1480-436-588

Mindray Medical France SARL
Europarc Creteil 1 allée des cerisiers 94000 CRETEIL, France
Voice: (33-1) 4513-9150
Fax: (33-1) 4513-9151

Mindray Medical Sweden AB
Rissneleden 136, 174 57 Sundbyberg, Sweden
Voice: (46) 8 555 54 121

Fax: (46) 8 555 54 101
Mindray Medical Rus Co.Ltd
2 Zvenigorodskaya Street, 13, bld 41, 123022 Moscow, Russia
Voice: (7-499) 553-6036
Fax: (7-499) 553-6039

Mindray Medical Colombia S.A.S
Room No.702, Torre UNIKA, Carrera 9 No.77-67, Bogota, Colombia
Voice: (57-1) 3130892 / 3210916

Mindray do Brasil – Comércio e Distribuição de Equipamentos Médicos Ltda.
Rua Tavares Bastos, 329, Perdizes, São Paulo - SP,
CEP 05012-020
Voice: (55-11) 3124-8026
Fax: (55-11) 3078-8035

Mindray Medical Mexico S. de R.L. de C.V.
Félix Parra # 175 Colonia. San José Insurgentes Sur Delegación Benito Juárez 03900 México, D.F.
Voice: (52-55) 5661-9450 / (52-55) 5662-6620
Fax: (52-55) 5662-6597

PT.Mindray Medical Indonesia
The East 11th Floor, Unit 08 Jl. Dr. Ide Anak Agung Gde Agung Kav. E 3.2 No.1 Jakarta 12950
Voice: (6221) 2902-7280
Fax: (6221) 2902-7283

Mindray Medical Egypt
Ground Floor-11 El Imam Metwally El Sharawi Street,
Sheraton Heliopolis - El Nozha, Cairo, Egypt
Voice: (2) 02 22691047
Fax: (2) 02 22691046

Mindray Medical Australia
ADD: 19 Harker St, Burwood VIC 3125, Australia
Voice: (61-3) 9830 6610
Fax: (61-2) 9550 2841

Mindray Medical Nigeria

ADD: 1004 Estate,Victoria Island, Lagos, Nigeria
Voice: (+234) 7030682053

Mindray Medical India Pvt. Ltd. (Delhi Office)
16th floor, Building No. 9-B, DLF Cyber City, DLF Phase- 3, Gurgaon- 122002, Haryana
Voice: (91-124) 4632-488
Fax: (91-124) 4632499


Mindray Medical India Pvt. Ltd. (Mumbai office)
B-404, City Point, Andheri Kurla Road, Andheri East, Mumbai, 400 059, India
Voice: (91-22) 4020-0000
Fax: (91-22) 4020-0011

Mindray Medikal Teknoloji Istanbul Ltd. Sti.
KÖYALTİ MEVKİİ 29 EKİM CAD.N.1 KUYUMCUKENT 2.NO.PLZ.3.KT. D.1- 9 34520 BAĞÇELİEVLER/YENİBOSNA, İSTANBUL, TÜRKİYE
Voice: (90-212) 482 08 77
Fax: (90-212) 482 08 78

Mindray Medikal Teknoloji Istanbul Ltd. Sti. Branch Office Ankara
TURAN GÜNEŞ BULVARI ALEKSANDER DÜBÇEK CAD. NO:18/A ÇANKAYA/ANKARA, TÜRKİYE
Voice: (90-312) 441 96 92
Fax: (90-312) 441 96 93

Mindray Medical Vietnam Co.,Ltd.
Unit 2, 4th Floor, "B" Building, Broadway Office Park, 102 Nguyen Luong Bang Blvd., Tan Phu Ward, Dist. 7, HCMC, Vietnam
Voice: (84) 8 5413-6970
Fax: (84) 8 5413-6971

Mindray Medical Thailand Limited
1768 Thai summit Tower Room 407, 4th floor, New petchaburee Road, Kwaeng Bangkok, Khet Huai Khwang, Bangkok 10320, Thailand
Voice: (66) 2652-7288
Fax: (66) 2652-7266

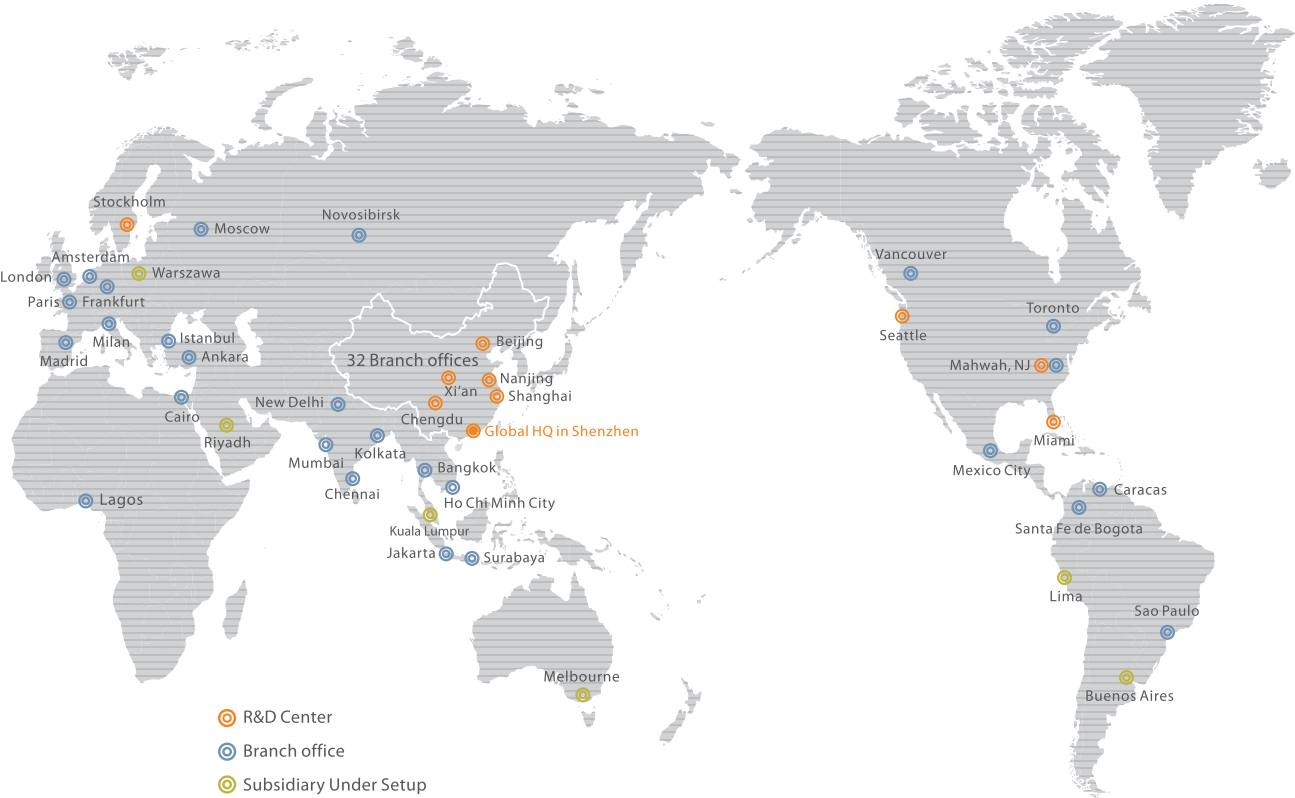


HyBase 8300/8500 Operating Table



mindray
healthcare within reach

<p>Mindray Building, Keji 12th Road South, High-tech Industrial Park, Nanshan, Shenzhen 518057, P.R. China Tel: +86 755 8188 8998 Fax: +86 755 26582680 E-mail: intl-market@mindray.com www.mindray.com</p> <p>Mindray is listed on the NYSE under the symbol"MR"</p>	<p>mindray healthcare within reach are registered trademarks or trademarks owned by Shenzhen Mindray Bio-medical Electronics Co., LTD. ©2014 Shenzhen Mindray Bio-Medical Electronics Co., Ltd. All rights reserved. Specifications subject to changes without prior notice. Pls refer to user manual for manufacturer information P/N:ENG-HyBase 8300/8500-210X285X16P-20150121</p> <div><div>CE</div><div>mindray</div></div>
--	--

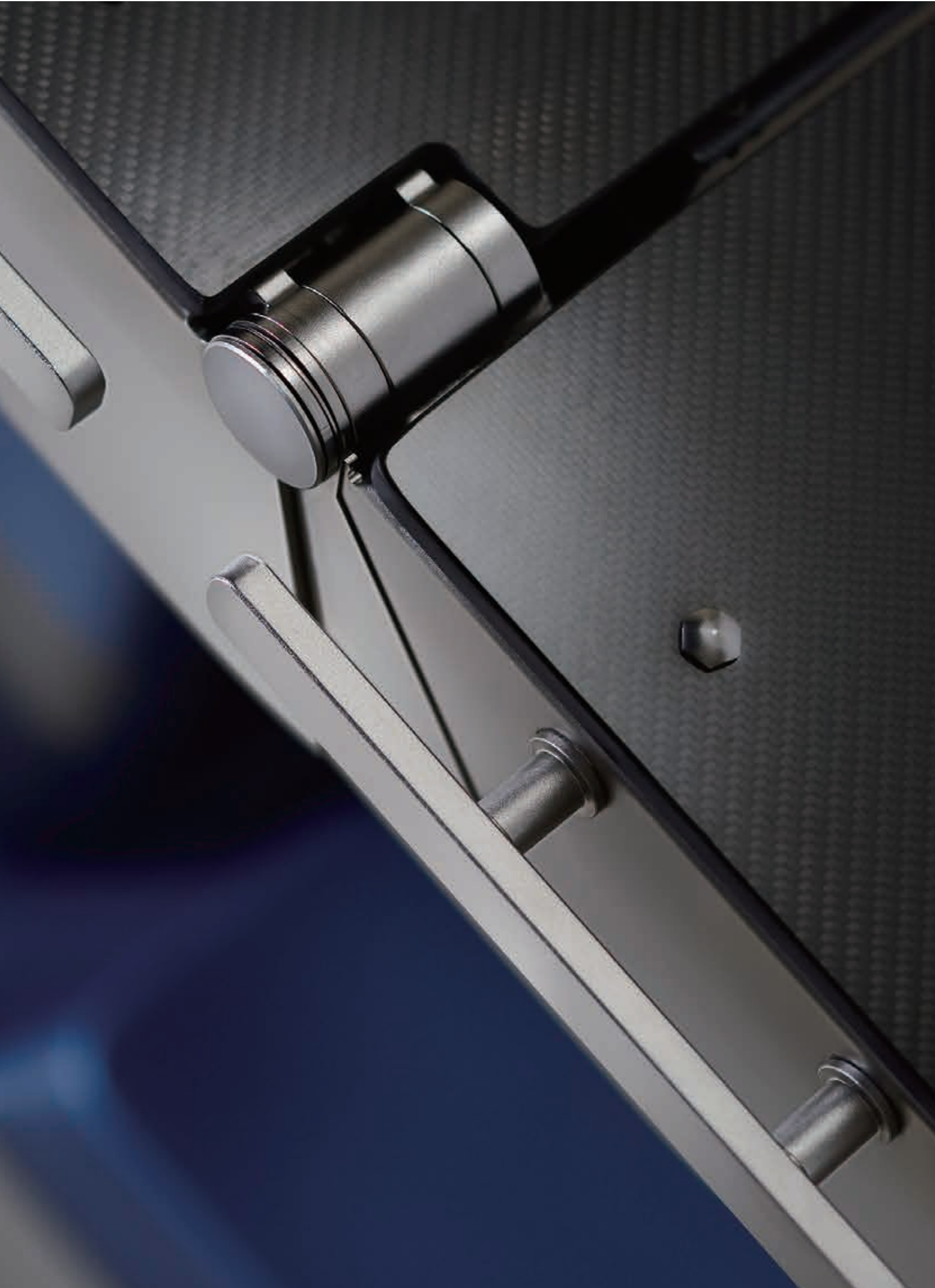

HyBase 8300/8500

Electro-Hydraulic Operating Table

Engineered for increasingly complex surgical scenes, HyBase 8300/8500 comes with extreme weight capacity up to 460 kg guaranteeing highest safety and stability in any position. With longitudinal slide up to 480 mm, it is ideally suitable for cardiovascular surgery with C-arm. The modular design tabletop enables it to be tailored for various surgical discipline needs.

- Longitudinal slide up to 480 mm for free access to C-arm
- Extreme weight load capacity of 460 kg with highest safety and stability
- Modular design tabletop tailored for various surgical discipline needs
- Intelligent anti-collision system to avoid table components collision during movement
- Multi-layer decompression pad with waterproof and seamless design
- 575 mm lowest table position for Neurosurgery
- IR wireless remote control system for free movements in OR

HyBase 8300/8500

Electro-Hydraulic Operating Table



Ergonomic and humanization design

- Soft material coating for comfortable holding
- Press any key to start and auto-power off
- Backlight for operating in a dark environment
- Flex, reflex, beach chair, and "0" position with one button for each position

01

Better visualization - 2.4" color LCD

The useful information will be shown on the display, such as movements, battery, brake, anti-collision information, etc.

02

03

10 memory table positions

The table can memorize 10 positions, so the user can select the memorized positions for the specific surgery freely without complex adjustment.

04

Intuitively UI design

The hand control button for table component movement is positioned in the corresponding patient figure chart, which enables the user to find the button fast and precisely.

05

IR wireless remote control system

With wireless hand control, the user can control the table in any place of the OR.

06

Wireless charging function

Without any charging interface, the hand control is easy to clean, matching the hygiene OR standard.

HyBase 8300/8500

Electro-Hydraulic Operating Table



Electric autodrives
casters for stable brake

Extreme Weight Capacity

HyBase 8300/8500 offers the highest level of safety and stability. With a maximum weight capacity of 460 kg, the table is designed to accommodate the lightest to the heaviest surgical patient in any surgery position, ideally suitable for bariatric patients.




Intelligent Anti-collision System

In the connection joint of leg plate, there is an intelligent module that can recognize each type of the table component and match the anti-collision data automatically. So the intelligent anti-collision system can monitor the positions of the table tops during movements such as sitting plate, leg plates and so on, to avoid table components collision.

Large Longitudinal Slide

With longitudinal slide up to 480 mm, HyBase 8300/8500 enables optimum access with the C-arm without repositioning the patient for imaging purpose, ideally suitable for such as cardiovascular surgery. By the way, the carbon-fibre table top (optional) provides excellent imaging quality.

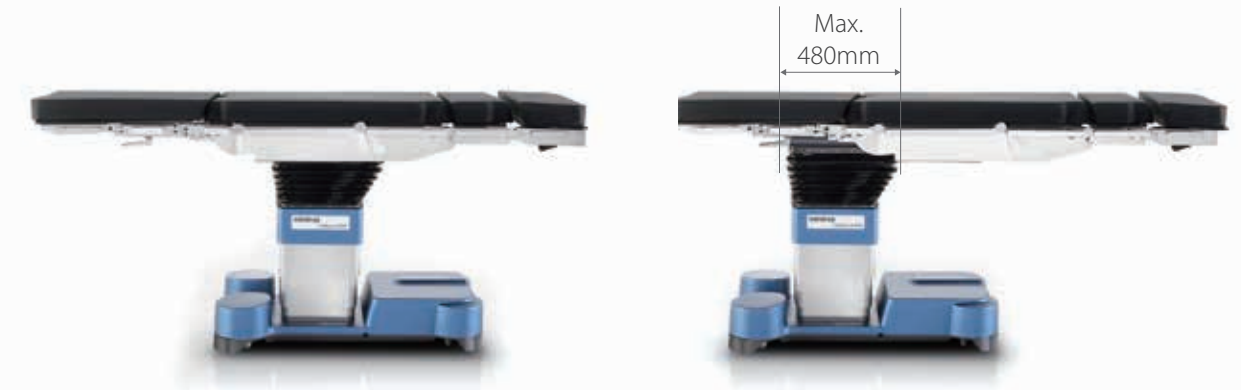

Multi-layer Decompression Pad

With the Special Foam Core pad and multi-layer design, HyBase 8300/8500 provides excellent pressure distribution, significantly avoids the main cause of decubitus, reduces the rate of postoperative health risk, and speeds up the patient recovery. Meanwhile, the pad is waterproof and easy to clean as ultrasonically-welded technology and GORE-TEX material.



Lowest Table Height

With lowest height up to 575mm for neurosurgery, the doctor can sit and do the surgery with comfort.




HyBase 8300/8500

Electro-Hydraulic Operating Table

Modular Design

Modular design table top is tailored for various surgical discipline needs, which also can low down the cost of the hospital as they only need to buy the specific module for more surgeries not one more table. With Module Recognition System, the table can recognize the module type and match movement data automatically.

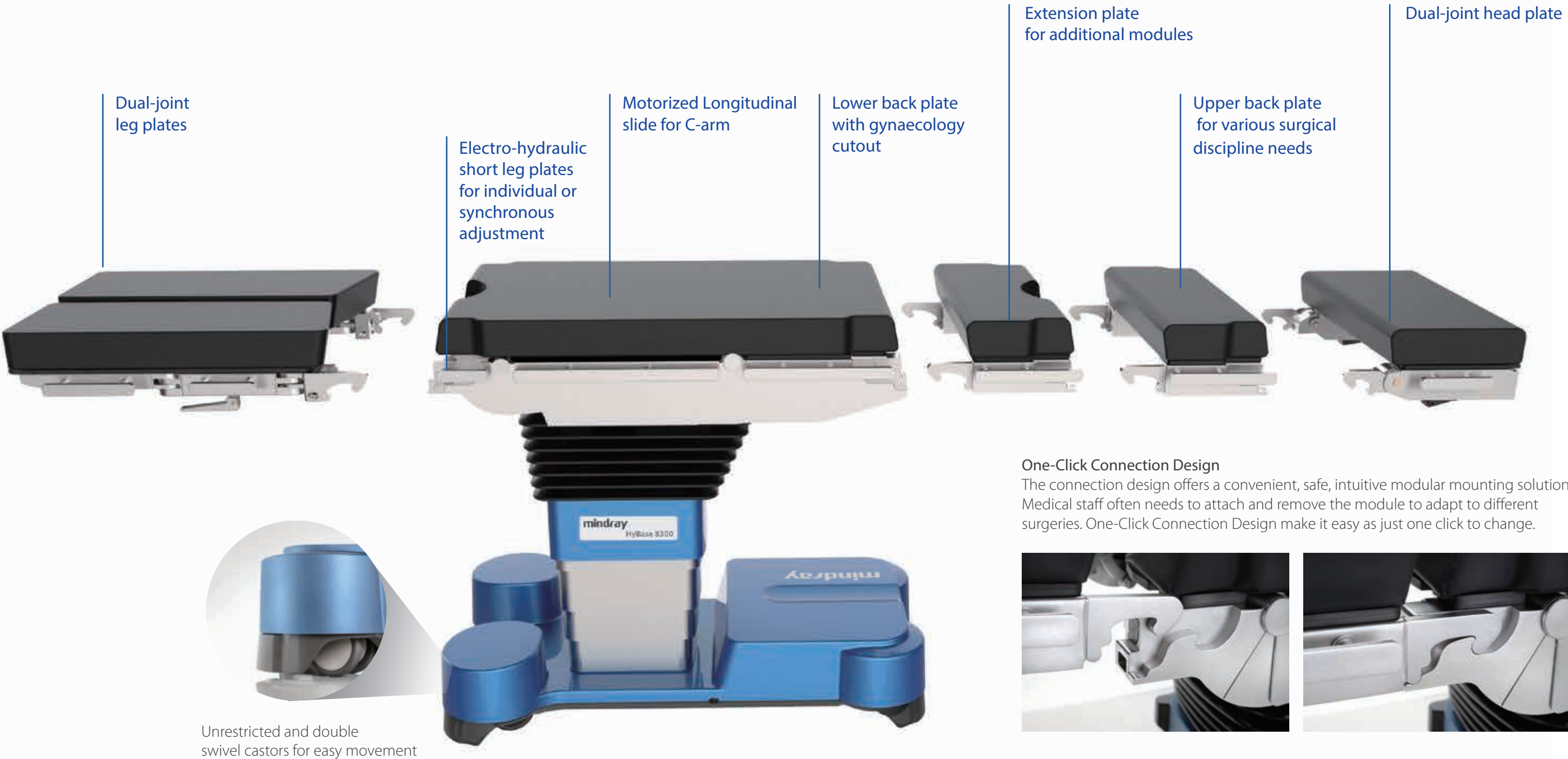

One-Click Connection Design
The connection design offers a convenient, safe, intuitive modular mounting solution. Medical staff often needs to attach and remove the module to adapt to different surgeries. One-Click Connection Design make it easy as just one click to change.





Cardiac OR




Orthopedics OR



Neurosurgery OR



Cardiovascular Surgery

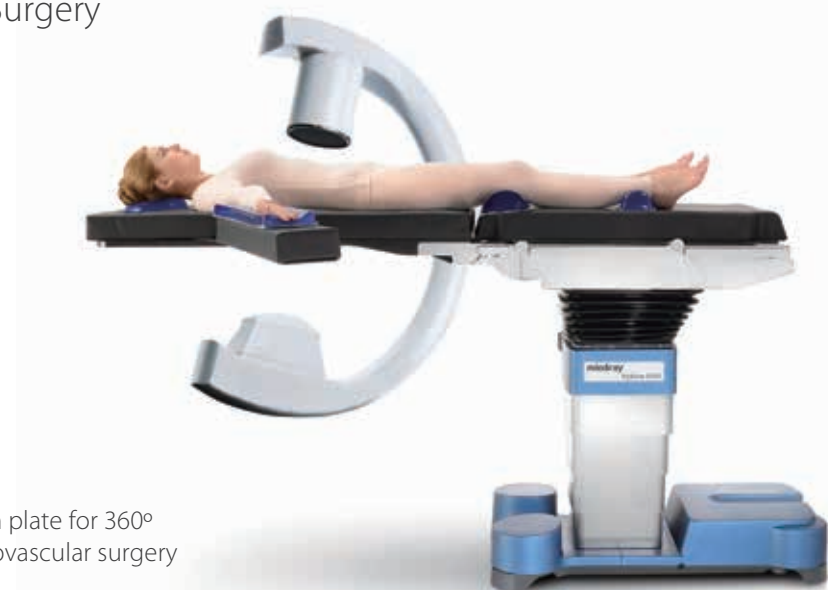


Carbon fiber extension plate for 360° access to C-arm, cardiovascular surgery

Orthopedic Surgery



Carbon fiber orthopedics extension for 360° access to C-arm

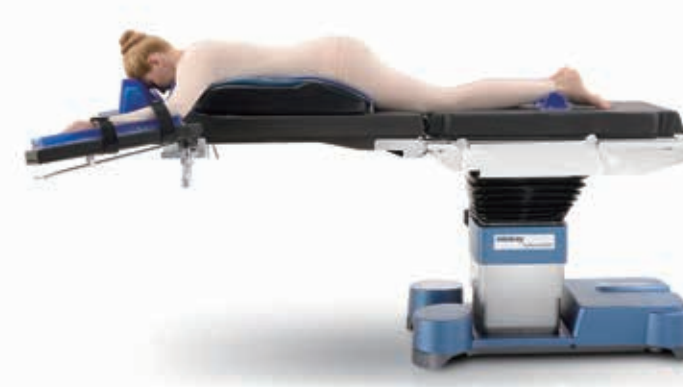

Carbon fiber extension plate for 360° access to C-arm, spinal surgery

Neurosurgery



Lowest table height with Skull Clamp (DORO) for Neurosurgery



Beach chair position for Neurosurgery

HyBase 8300/8500

Electro-Hydraulic Operating Table

Flexible Patient Positioning

General Surgery



Struma Surgery



Gynecology and Urology Surgery



Reverse mode, after sliding to head end



Reverse mode, with extension plate



Reverse mode, ergonomic positioning with stirrups

ENT Surgery



Lowest table height for ENT surgery

Kidney Surgery



Lateral position with dual-head plate

Endoscopic Surgery




More access to patient with dual-joint leg plate

HyBase 8300/8500


Electro-Hydraulic Operating Table

Technical Specifications


Dimensions	
Length (Main unit)	970mm
Height without padding	HyBase 8300: 575-1,075mm/HyBase 8500: 700-1,200mm
Width without side rail	520mm
Width with side rail	590mm

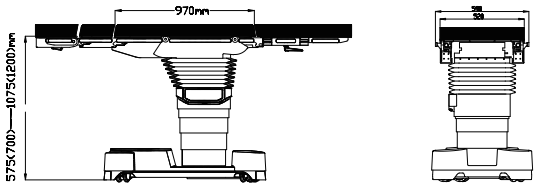


Trendelenburg position	
Trendelenburg	HyBase 8300: Max. 36° / HyBase 8500: Max. 46°
Reverse trendelenburg	HyBase 8300: Max. 36° / HyBase 8500: Max. 46°

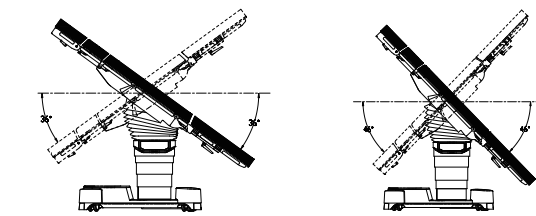


Leg plates	
Leg plates up	Max. 80°
Leg plates down (Without sliding)	Max. 90°
Leg plates down (With sliding)	Max. 100°

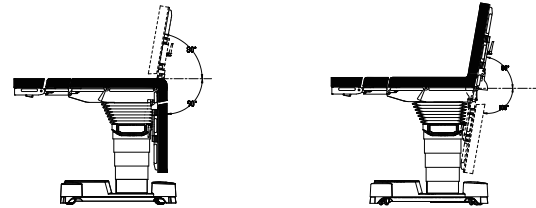


Reflex Position	
Reflex Position	110°
Trendelenburg	30°
Mounting point up	40°

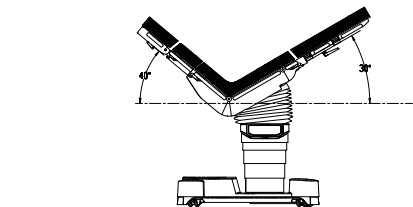


Longitudinal slide (optional)	
Longitudinal slide	Max. 480mm

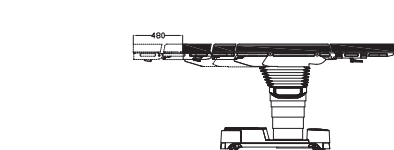


Lateral tilt	
Left tilt	HyBase 8300: Max. 26° / HyBase 8500: Max. 36°
Right tilt	HyBase 8300: Max. 26° / HyBase 8500: Max. 36°

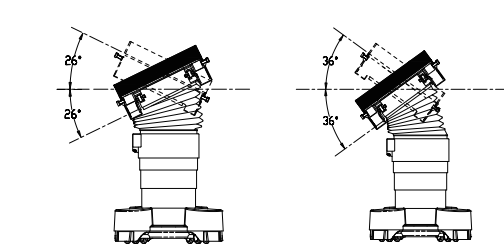


Back plate	
Back plate up	Max. 90°
Back plate down	Max. 45°

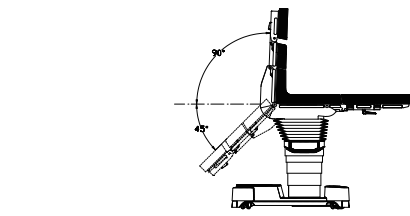

Head plate	
Head plate up	Dual-joint: Max. 60°/Single-joint: Max. 45°
Head plate down	Dual-joint: Max. 60°/Single-joint: Max. 90°

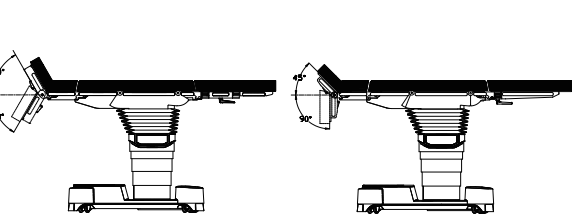

Flex Position	
Flex Position	220°
Reverse Trendelenburg	20°
Mounting point down	20°

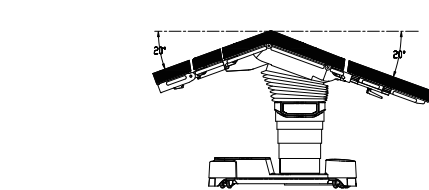

Beach-chair position	
Back plate up	80°
Trendelenburg	25°
Leg plates down	70°

