

## DC-60 Pro

Diagnostic Ultrasound System

An all-around partner to trust

X-Insight



# X-Insight

## An insightful solution to envision more.

Mindray's brand-new solution is an excellent fusion of continuous customer insights into clinical needs, combined with continually evolving cutting-edge ultrasound technologies.

Full of vitality, with the goal of insight into the future and pursuit of the infinity, it is constantly improving with increased scalability.

## High efficiency with precision imaging

As an all-around partner, DC-60 Pro with X-Insight devotes to provide a comprehensive solution to help you manage all-aspect daily clinical practice with ease and certainty.

Based on the deep insights of customer needs, the DC-60 Pro with X-Insight is designed to deliver high efficiency with precision imaging, which is empowered by eXpress Clarity, eXceptional Intelligence and eXceeding Experience :



eXpress Clarity



eXceptional Intelligence



eXceeding Experience



# eXpress Clarity


More clarity at hand

To achieve excellent images with minimum effort is your highest priority and our endless pursuit. With continuous innovation in imaging and transducer technologies, the DC-60 Pro with X-Insight delivers immediate clarity with ‘Touch and Bingo’ to allow you to get the optimum images that you need immediately, as soon as the transducer touches the body.

## ‘Touch and Bingo’ with evolutionary transducers


### Optimized transducers with 3T technology

Combining with Mindray unique 3T technology (Triple-matching layers, Total-cut design, Thermal control), the optimized convex, phased array and volume transducers provide a wider range application coverage, resulting in an optimum scanning solution in ABD, Cardiology, OB/GYN and more.

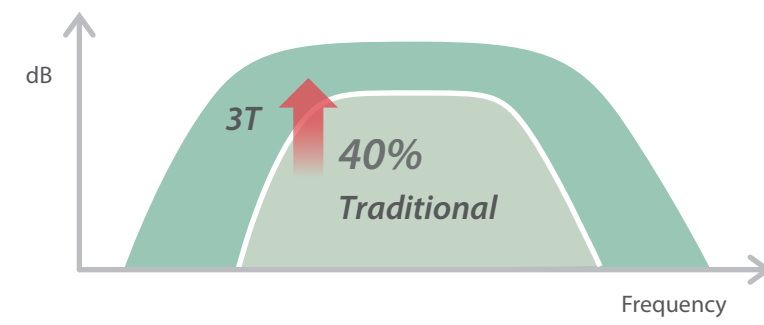

### ComboWave transducer with 3T technology

Compared with traditional transducers, ComboWave transducers utilize a new type of composite piezoelectric material to dramatically optimize the acoustic spectrum and reduce acoustic impedance. Further integrated with Mindray unique 3T technology, the ComboWave linear transducers allow you to experience outstanding performance with extreme image resolution and uniformity in thyroid, breast, vascular, and more.

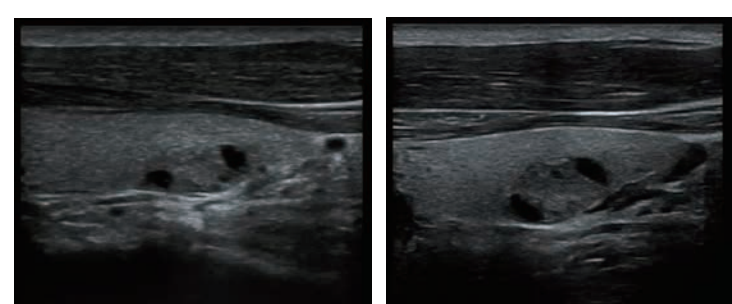


Traditional Linear Image

ComboWave Linear Image

## Thunder-speed imaging powered by X-Engine

The new X-Engine utilizes the evolutionary integration of GPU and CPU to enable multi-core parallel processing for fast imaging. With the advanced imaging engine, the imaging processing speed is accelerated three or four times, resulting in extremely fast imaging for 3D/4D and other applications.

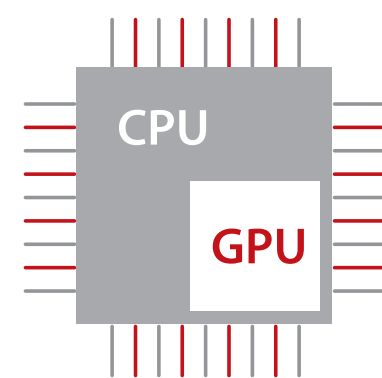

Processing Speed 3X Plus


## Ultra-Speed data transmission with security guarantee

By the optional integration of Solid-State Disk(SSD), DC-60 Pro with X-Insight provides the ultra-speed data transmission capability to accelerate the system operational speed, for instance, faster imaging, quicker boot-up or wake-up. And the SSD is able to extend the data life-time of ultrasound system, which enhances the protection of patient data and privacy guarantee.




Transmission speed 5X Plus

## Proven advanced applications:

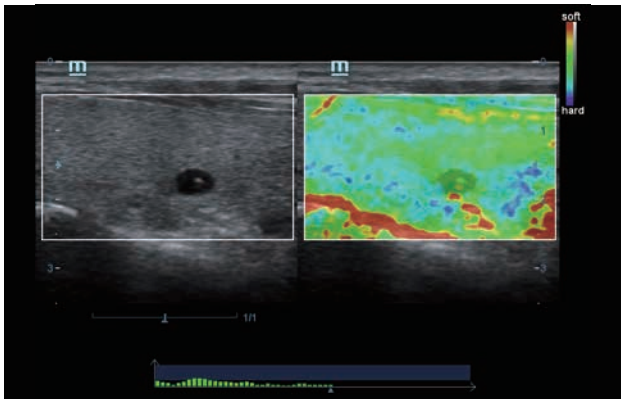

Natural Touch Elastography

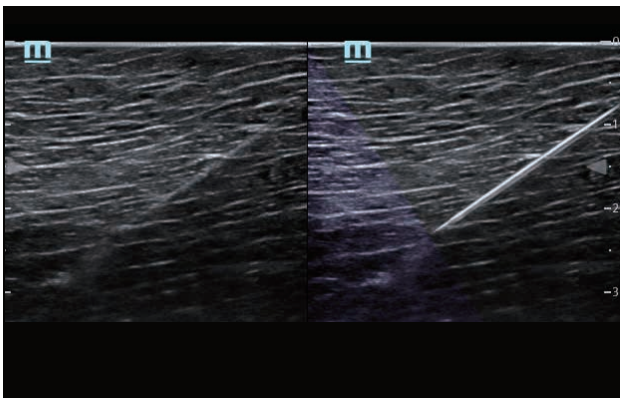

iNeedle+ off

iNeedle+ on




# Image gallery




Liver, 2D



Kidney HR Flow

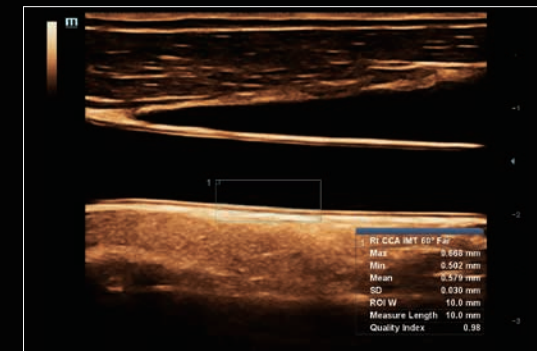


Common Carotid Artery IMT

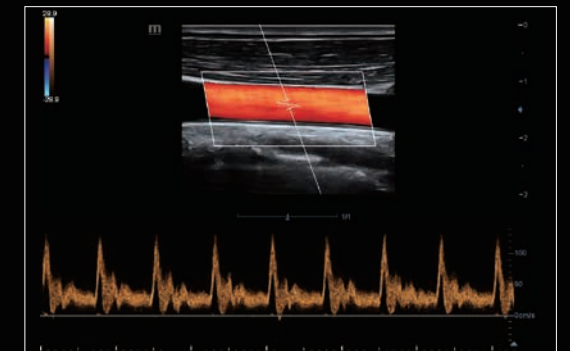


Common Carotid Artery PW




Abdomen Mass

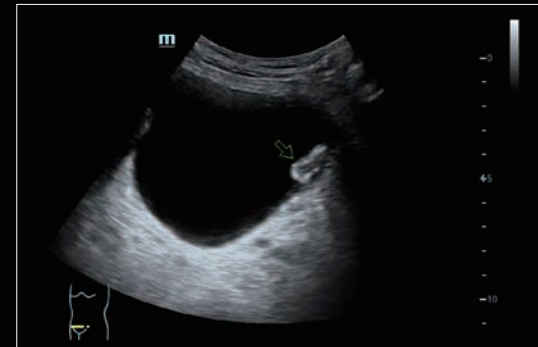


Bladder Mass, 2D

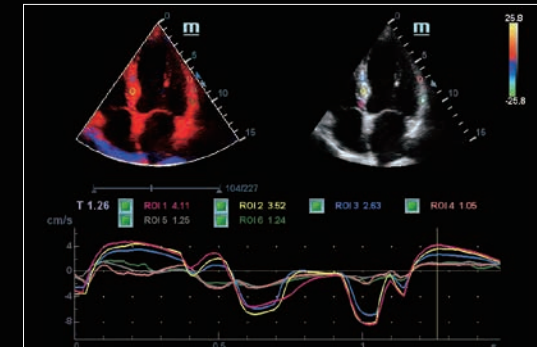


Cardiac TDI QA

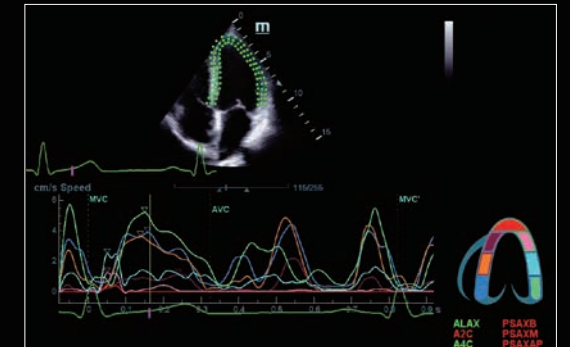


Cardiac TT QA

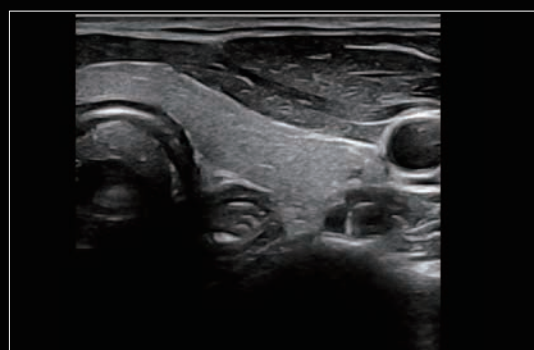

Thyroid, 2D

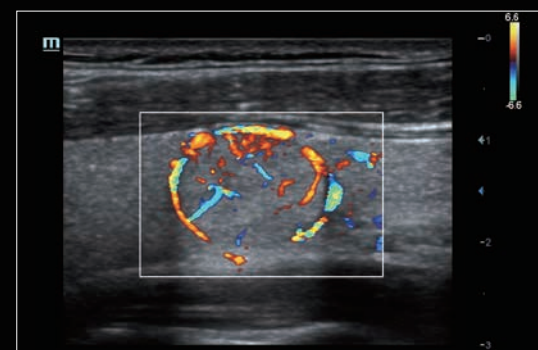


Thyroid Mass CFM

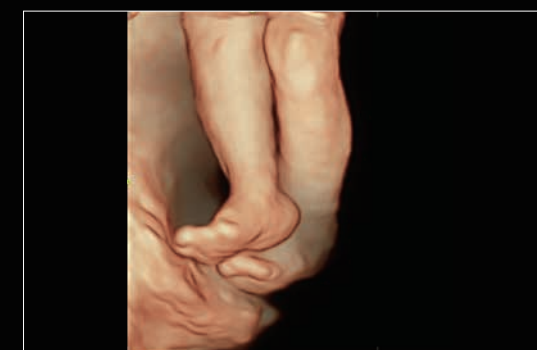


Fetal Leg 3D



Fetal Face iLive

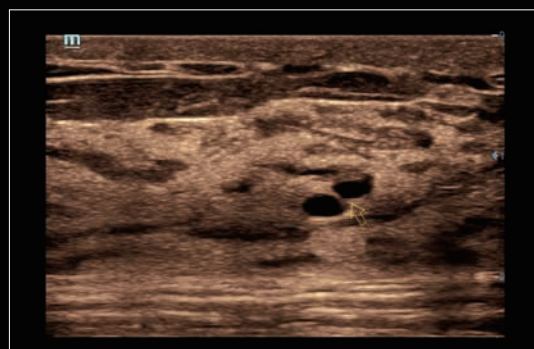


Dilation of Mammary Duct

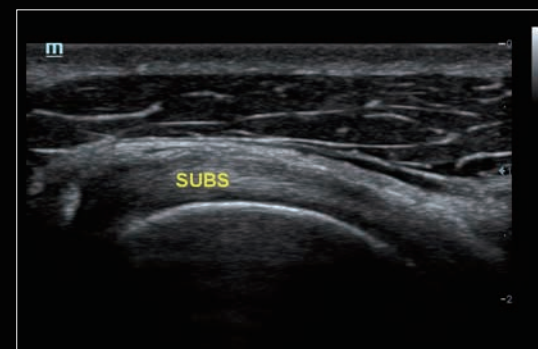

Subscapular Muscle Tendon



Fetal Heart iLive

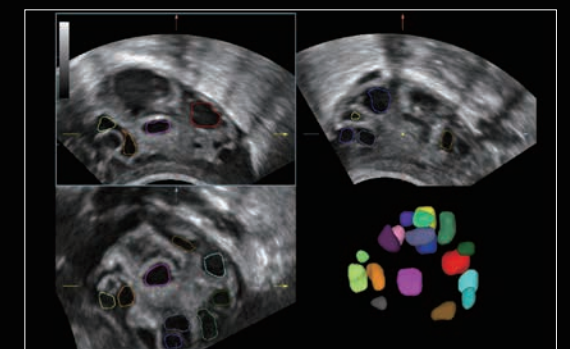

Smart Follicle






## Smart Planes CNS

More accuracy with smart acquisition

Smart Planes CNS is a user-friendly solution to improve through-put, and reduce user dependency. **With a simple button-click** on a 3D fetal brain volume image, the standard CNS scanning planes (MSP, TCP, TTP and TVP) and a range of related anatomical measurements (BPD, HC, OFD, TCD, CM and LVW) are obtained immediately with high accuracy, which reduces about 90 Second per exam.

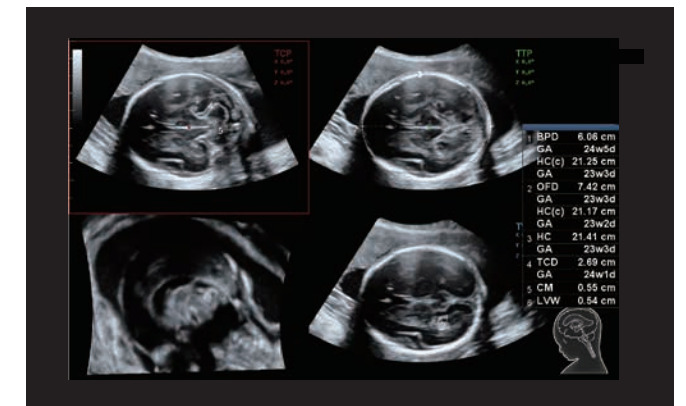


Smart Planes CNS

## Smart Track

More simplicity with auto image optimization

Smart Track provides fast and intelligent optimization for vascular imaging with simply **one-touch operation**. It can optimize Color, Power and PW spectrum by **self-tracking** and reduce time-consuming steps. Therefore the vascular exam workflow is simplified with an optimal view.

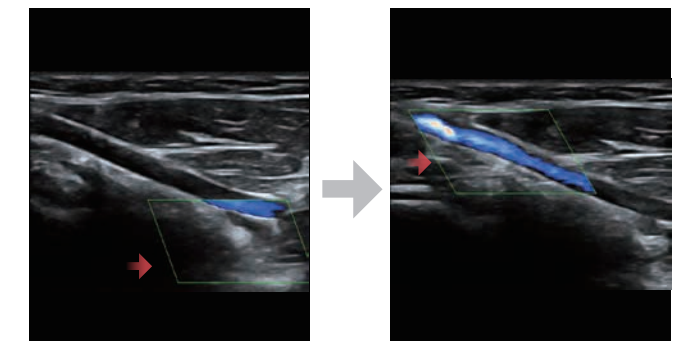


Smart Track

## iWorks

More consistency with exam protocol

iWorks takes advantage of the inbuilt standard scan protocols for more consistency and **reduces the exam time by up to 50%**. The highly flexible and powerful feature is further enhanced with user-defined capabilities.

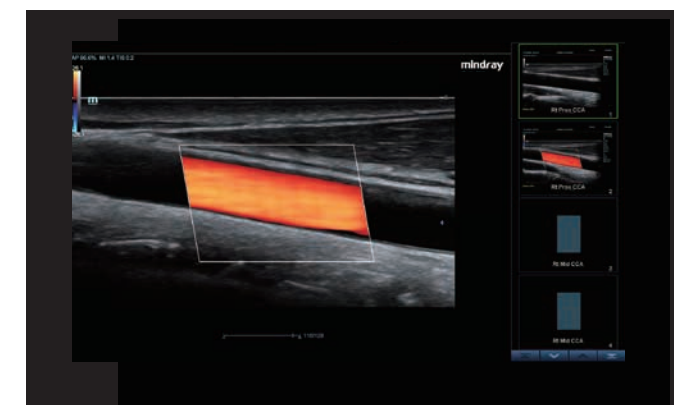

iWorks


# eXceptional Intelligence

Intelligence throughout entire workflow

To improve scanning efficiency with more accuracy and consistency, the DC-60 Pro with X-Insight provides All-Smart Exams for all kinds of advanced applications, from plane acquisition to image optimization, and from calculation to exam protocol.

## Smart Face

Smart Face provides fast and intelligent optimization for fetal face with one-touch operation. It immediately removes occlusions and eliminates noise information, generating an optimal view of the fetal face with more simplicity.

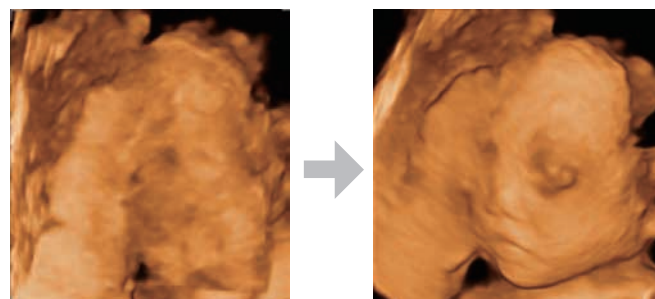


Smart Face

# eXceeding Experience

## Experience with high productivity

Comfortable experience in scanning results in more patient focus for high quality exams. The DC-60 Pro with X-Insight provides outstanding ease of use with better ergonomics, easier exams, and flexible management, even beyond your expectations.

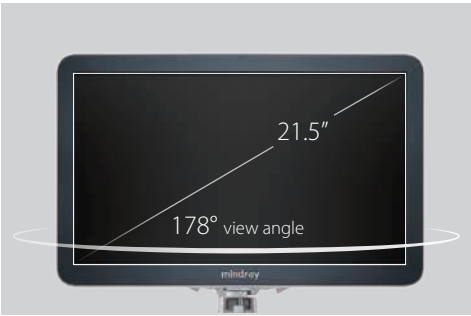

21.5" full HD monitor



Unique dual-wing floating arm(optional)



13.3" ultra-slim touch screen with angle adjustment



Multi-Gesture operation



Interactive ultrasound APP  
--- MedSight/MedTouch



Elevating and rotating control panel

## Higher Reliable Dependence

For decades, Mindray cares about using intelligence to enhance the safety, accuracy and stability of its products, to ensure the satisfaction during your daily diagnosis.

### Guarantee of safety and stability:



**EMC (Electro Magnetic Compatibility)**  
By obtaining the CB Test Certificate (the certificate for EMC approval, Ref. Certif. No. FI-31332), DC-60 Pro with X-Insight reaches the highest anti-interference level assurance, which ensures the stability of its imaging quality.




**Class B (Power supply requirement)**  
The power supply level required by DC-60 Pro with X-Insight is Class B. It can work under not only hospital power supply condition (Class A), but city power supply condition (Class B), and it offers you more safety and compatibility when you operate DC-60 Pro with X-Insight in any complicated environment.



Integrated control panel,  
anti-dust design



Transducer cable management,  
designed for regular and protection

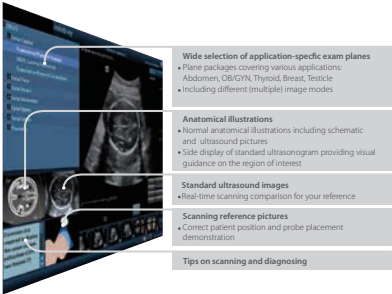

iScanHelper (built-in tutorial)

# May the X-Insight be with you!

As an all-around partner to trust, DC-60 Pro with X-Insight will provide you 'high efficiency with precision imaging' from all-aspects of your daily clinic application.